


## II. Installation and preparation

1. Check product's composition
2. Check before product installation
3. Installation
4. Check after installation



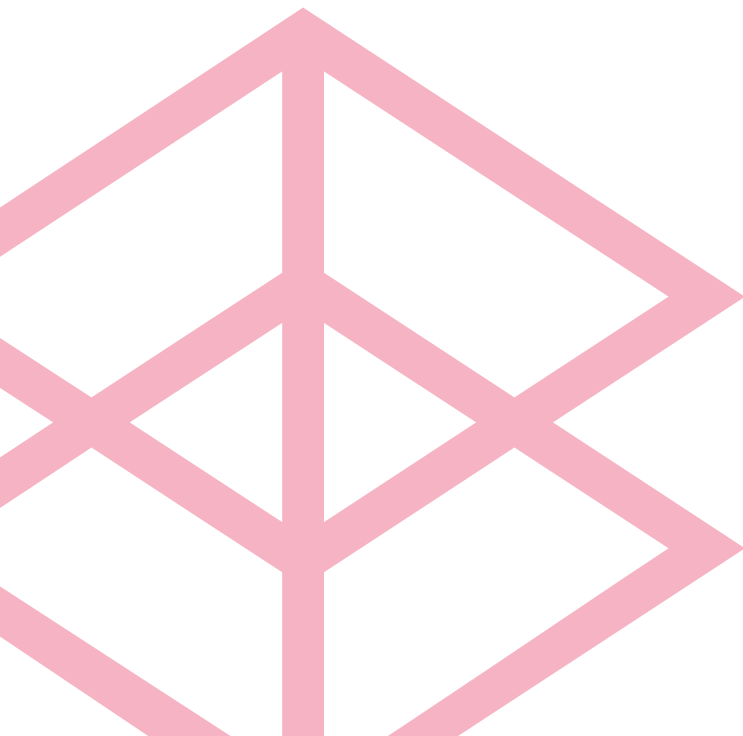

# 1. Check product's composition

1		<b>External projectors (3-Sphere, 9-Sphere)</b> Easy to use on body parts such as the abdomen, shoulders, arms, legs and etc. Can be used with the main product or shared between two people.
2		<b>External Projector Cable</b> To use the external projector, connect it to the main body.
3		<b>External projector cover</b> It can be used while covering the external projector.
4		<b>Remote control</b> Used to check the product functions and operating condition of the product.
5		<b>Remote Control Cable</b> To use the remote control, connect it to the main body.
6		<b>Remote control cradle</b> The remote control cradle can be assembled onto and used on the main body of the product.
7		<b>Remote Control Back Label</b> Indicates the functions of the keys of the remote control; it can be attached to the back of the remote control.
8		<b>Power cord</b> Connects with the main body and plugs into an outlet for use.
9		<b>User manual</b> Information on how to use the product.
10		<b>SD card</b> Users may add their preferred music in addition to the Ceragem default sound.
11		<b>Outer fabric</b> Protects against product contamination and burns; this detachable fabric makes washing easy.
12		<b>Head cushion</b> This is attached to the top of the product.
13		<b>Assembly guide</b> Product assembly methods are described.



## 2. Check before product installation

1. Install where the ground is flat to avoid the risk of damaging the product.

2. Try and avoid installation where there is a lot of humidity or dust.  
(To avoid risking electric shocks or fire.)

3. Install according to the designated order.

4. When installing, check that the rated voltage is 220-240V~.  
(To avoid circuit damage or fire.)

5. Do not drop the product or apply any force to it.  
(This could cause malfunctioning.)

6. Do not use damaged power cords, plugs or old sockets.


7. You must connect to a socket with a connection terminal.  
(To avoid risking electric shocks.)



# 3. Installation


1

Take out the product and its components from the box.

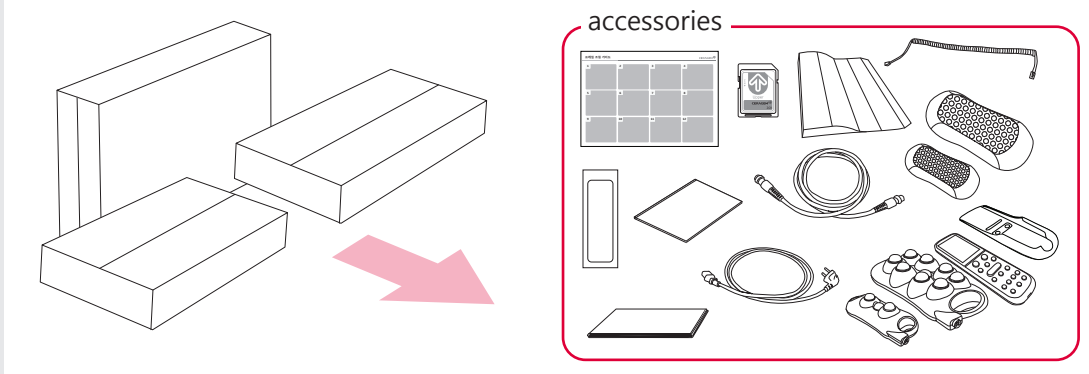

2


After fastening the supporting mat to the sliding rail, place the main body mat (which is attached to the head cushion and the outer fabric) onto the main body frame. (After moving the supporting frame to the middle of the main body, push in the main body mat completely.)

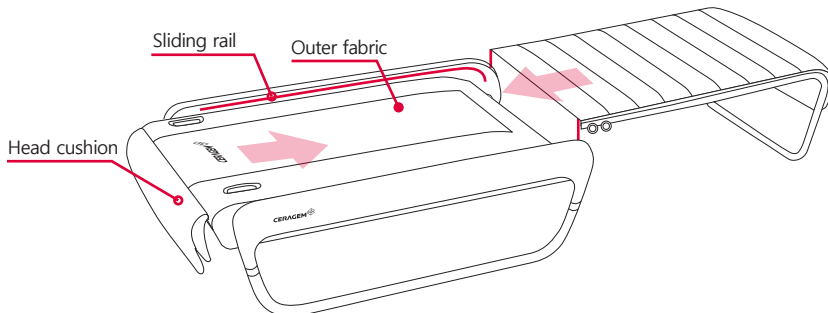

\* Please refer to the Frame Assembly Guide on how to assemble the frame.  
\* Please check both sides of the rail carefully before assembling the supporting mat.

3

Connect the cable which is attached to the main body to the connector of the supporting frame. (Open the supporting frame before connecting.)

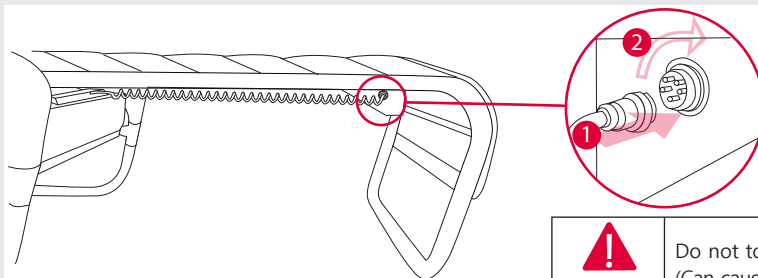


Do not touch the connector pin.  
(Can cause bodily harm or mechanical failure)




## 3. Installation

- 4 Please connect the Remote Control, External Projectors and cable properly.

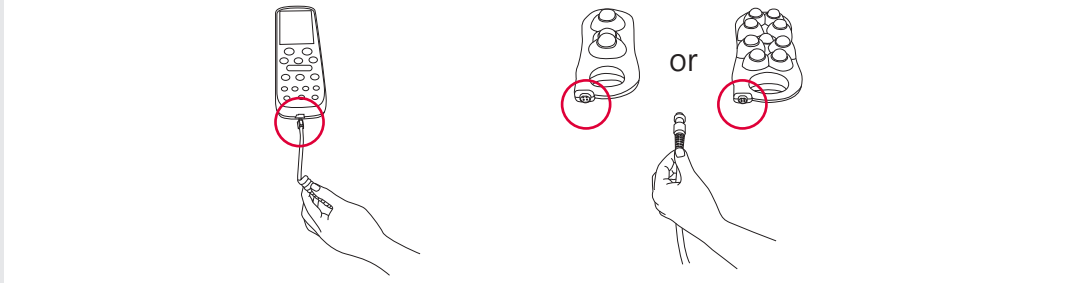

- 5 Assemble the remote control cradle onto the frame and place the remote control onto the cradle. Check the connectors for the SD card, remote control and the external projector and connect appropriately.

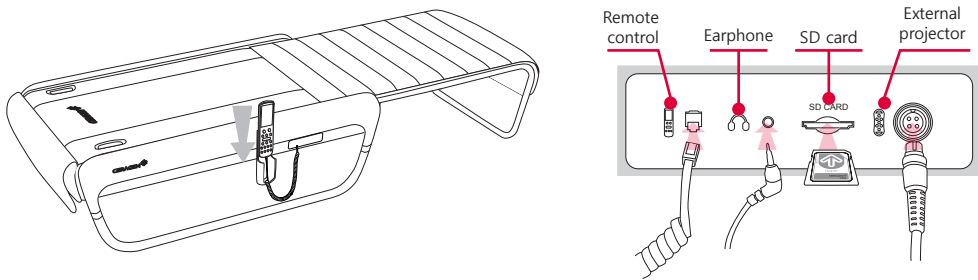

Caution


Please keep earphones plugged in only when in use.  
Please make sure that the sticker side of the SD card is facing up before inserting into the main body.  
After turning the power off, please push the SD card all the way in.



Warning


Do not touch signal input, signal output or other connectors, and the patient simultaneously.  
Do not touch the connector pin.(Can cause bodily harm or mechanical failure.)

- 6 Connect the power cord onto the main body and plug into a socket..

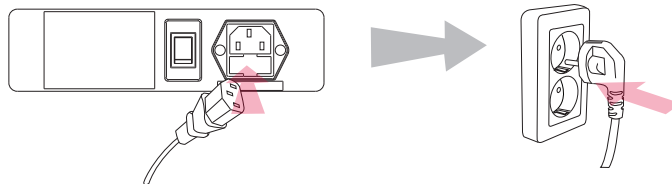

Caution


If the power cord is connected when the power switch is "On", it can cause mechanical failure.

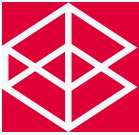

## 4. Check after installation


1

Turn on the power at the side of the main body.  
(Please turn off the power when the product is not being used.)

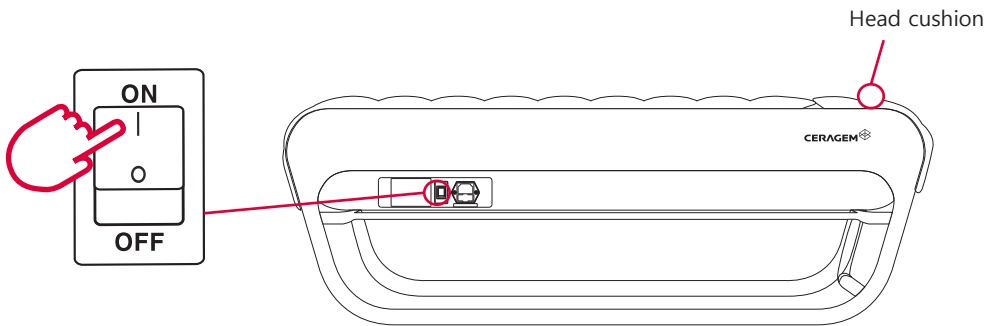

2

Press the power key of the remote control and check whether the display turns on.

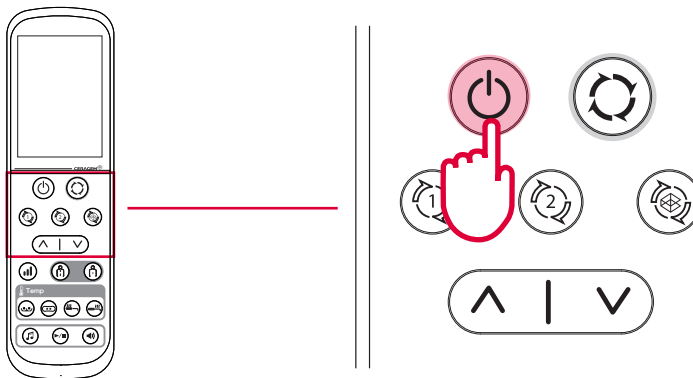

\* For the remote control's screen, refer to page 20.



Caution

1. Before connecting the power cord, connect the remote control first.
2. According to the user, set the temperature appropriately when using the product.
3. When sitting down and using the product, open the supporting frame completely.
4. After using the product, turn off the power on the remote control and then turn off the main body.